

TECHNICAL MANUAL

INSTRUCTIONS FOR PROCUREMENT, ISSUE, USE AND MAINTENANCE OF MEDICAL KITS

Prepared By: Automated Technical Order System (ATOS)

THIS PUBLICATION SUPERSEDES TO 00-35A-39, DATED 15 JUNE 2000. THIS PUBLICATION HAS BEEN EXTENSIVELY REWRITTEN, THEREFORE REVISION BARS HAVE BEEN OMITTED.

DISTRIBUTION STATEMENT - Approved for public release; distribution is unlimited.

LIST OF EFFECTIVE PAGES

NOTE: The portion of the text affected by the changes is indicated by a vertical line in the outer margins of the page. Changes to illustrations are indicated by miniature pointing hands. Changes to wiring diagrams are indicated by shaded areas.

Dates of issue for original and changed pages are:

Original 0..... 1 April 2001

TOTAL NUMBER OF PAGES IN THIS PUBLICATION IS 42, CONSISTING OF THE FOLLOWING:

Page No.	*Change No.	Page No.	*Change No.	Page No.	*Change No.
Title	0				
A.....	0				
i - ii	0				
1 - 38.....	0				

*Zero in this column indicates an original page

TABLE OF CONTENTS

Section/Para	Page
1 Purpose.....	1
2 Authorization for Use and Issue Procedures.....	1
3 Controlled Substances.....	1
4 Marking of First Aid Kits.....	2
5 Kit Components.....	2
6 Serviceability Tags.....	3
7 Inspection of Medical Kits.....	3
8 Inspection By Other Than Medical Personnel.....	4
9 Funding for Replacement Components.....	4
10 Turn In of Excess Kits.....	4
11 Furnishing Kits to Air Force Contractors.....	5
12 Kit Standardization or Modification Procedures.....	5
13 NSN 6545-00-116-1410, First Aid Kit, General Purpose, Rigid Case.....	5
14 NSN 6545-00-139-3671, Survival Kit, Individual.....	5
15 NSN 6545-00-782-6412, Survival Kit, Individual, Tropical, Tactical Aircraft.....	14
16 NSN 6545-00-853-6309, First Aid Kit, Eye Dressing.....	20
17 NSN 6545-00-919-6650, First Aid Kit General Purpose, Flexible Case, Panel Mounted.....	20
18 NSN 6545-00-922-1200, First Aid Kit General Purpose, Rigid Case.....	24
19 NSN 6545-00-927-3000, Medical Instrument and Supply Set, Aeronautic Emergency.....	27
20 NSN 6545-01-090-0645, First Aid Set, Shelter.....	27
21 NSN 6545-01-094-8412, First Aid Kit, Individual.....	30
22 NSN 6545-01-400-3397, First Aid Kit, Individual, Modified.....	31
23 NSN 6545-01-120-2632, Survival Kit, Individual, Tropical, Tactical Aircraft, in a Plastic Bag.....	34
24 Component Cross Reference to First Aid Kits.....	36

LIST OF ILLUSTRATIONS

Figure	Title	Page
1	AFTO Form 104 and DD Form 1574.....	3
2	First Aid Kit, General Purpose, Rigid Case, NSN 6545-00-116-1410	6
3	Kit Assembly.....	16
4	Survival Kit, Individual, Tropical, Tactical Aircraft, NSN 6545-00-782-6412.....	19
5	First Aid Kit, General Purpose, Flexible Case, Panel Mounted, NSN 6545-00-919-6650	23
6	First Aid Kit, General Purpose, Rigid Case, NSN 6545-00-922-1200	26
7	First Aid Kit, Individual, NSN 6545-01-094-8412	33

LIST OF TABLES

Number	Title	Page
1	Special Kits	1
2	NSN 6545-00-116-1410 - Components	8
3	NSN 6545-00-139-3671 - Components	11
4	NSN 6545-00-782-6412 - Components	17
5	NSN 6545-00-853-6309 - Components	20
6	NSN 6545-00-919-6650 - Components	21
7	NSN 6545-00-922-1200 - Components	25
8	NSN 6545-00-927-3000 - Components	28
9	NSN 6545-01-090-0645 - Components	29
10	NSN 6545-01-094-8412 - Components	31
11	NSN 6545-01-400-3397 - Components First Aid Kit (Modified Kit of NSN 6545-01-094-8412).....	32
12	NSN 6545-01-120-2632 - Components	35
13	Component Cross Reference to First Aid Kits	37

1 PURPOSE.

- 1.1 To list special purpose medical first aid and/or survival kits authorized for use by USAF activities.
- 1.2 To describe the procedures to standardize or modify medical first aid or survival kits.
- 1.3 To delineate responsibilities for authorizing, requisitioning, issuing, installing, inspecting, controlling and disposing of medical kits.

2 AUTHORIZATION FOR USE AND ISSUE PROCEDURES.

2.1 Authorization for issue of special kits. (See table 1.)

2.2 In addition to the authorizations listed below, various USAF technical orders (TOs) list one or more of these kits for installation in aircraft, survival kits, or other major items. Organizations requiring issuance of medical kits not specifically authorized by a current Table of Allowance (TA) will cite the appropriate TO as the authority for issue.

2.3 Special Allowance Source Code (ASC) 041 authorizations may be granted by the Director of Base Medical Services (DBMS) as specified in AFMAN 23-110, Volume 5, chapter 10, paragraph 10.4.12.

2.4 Items listed in this TO will be initially issued to using activities on a reimbursable basis by the local Medical Materiel Manager (MMM). Requests for initial issues will be submitted to the MMM through the local Equipment Management Office (EMO) for verification of TO or TA authorization. Requests for replacement kits may be forwarded directly to the MMM. Initial issue of kits to USAF contractors is discussed in paragraph 11.

2.5 The kits listed are designed for the specific purpose and use intended and will not be modified, disassembled, issued, or used in any manner not prescribed herein, which includes for practice or play exercises. However, if required to keep flight crews mission capable, the Medical Group Commander may approve temporary substitutes after consulting with the senior flight surgeon in the Wing, or on the base, preferably an Aerospace Medicine Specialist. Kits containing these will be marked with an expiration date of 30 days, during which time permanent substitutes will be coordinated with AFMLO/FOC-D and the Air Force consultant for first aid kits. Medical kits will not be used when other medical services and supplies are readily available. The kits need to be ready and completely assembled for their intended purpose.

3 CONTROLLED SUBSTANCES.

All narcotics and non-narcotic controlled substances have been deleted from these medical kits. Controlled substances are normally not included in medical kits. Amphetamines will not be added to medical kits under any circumstances.

Table 1. Special Kits

NATIONAL STOCK NUMBER	ITEM DESCRIPTION	ASC/REMARKS
6545-00-116-1410	First Aid Kit, General Purpose, Rigid Case	016B, 410, 429, 456, 538
6545-00-139-3671	Survival Kit, Individual	041
6545-00-782-6412	Survival Kit, Individual, Tropical, Tactical Aircrew	041
6545-00-853-6309	First Aid Kit, Eye Dressing	Component part of larger assemblies, see component parts list elsewhere in this TO
6545-00-919-6650	First Aid Kit, General Purpose, Flexible Case, Panel Mounted	016B

Table 1. Special Kits - Continued

NATIONAL STOCK NUMBER	ITEM DESCRIPTION	ASC/REMARKS
6545-00-922-1200	First Aid Kit, General Purpose, Rigid Case	016B, 156, 538, 699
6545-00-927-3000	Medical Instrument and Supply Set, Aeronautic Emergency	016B
6545-01-090-0645 6545-01-094-8412	First Aid Set, Shelter First Aid Kit, Individual: Packed in Case, 6545-01-094-6142	Buddy Care and/or Self Aid Supplies 016B, 538
6545-01-120-2632	Survival Kit, Individual, Tropical, Tactical Aircrew, In Plastic Bag	041
6545-01-400-3397	First Aid Kit, Individual (modified kit of 6545-01-094-8412)	016B, 538

4 MARKING OF FIRST AID KITS.

4.1 To permit rapid identification of first aid kits by crew members and inspectors, the following kits will be marked as prescribed in paragraph 4.1.1 prior to installation in an aircraft, vehicle, or water craft.

6545-00-116-1410	6545-00-922-1200
6545-00-919-6650	6545-00-927-3000
6545-01-090-0645	

4.1.1 One symmetrical red cross on a white square background not less than 1 nor more than 4 inches square (depending on available space) will be painted or otherwise permanently applied to the side or end of each kit where it is readily visible.

4.1.2 If a kit is installed in a container bracket or holder which does not permit visibility of applied markings, an additional marking will be applied by base maintenance personnel. This marking will be placed on the container, bracket or holder, or on a structural part adjacent to the kit.

4.2 Special materials suggested for use in marking are listed below:

4.2.1 Red Cross Symbol 3 X 3 inches on 4 X 4 inches white background, WR-ALC, part number 3613053. Requests for this item are to be submitted on AF FORM 764A, REQUISITION AND REQUIREMENT REQUEST, to the 2853rd Air Base Group, DAPT, Robins AFB, GA 31098.

4.2.2 NSN 8010-00-141-2952 Lacquer, Nitro-cellulose, Red (Source GSA).

4.2.3 NSN 8010-00-290-6983 Lacquer, Nitro-cellulose, White (Source GSA).

5 KIT COMPONENTS.

5.1 The components of each kit listed in this T.O. are listed in the Universal Data Repository (UDR). In the event discrepancies exist between this T.O. and the catalog, using activities are to immediately notify AFMLO/FOC-D, 1423 Sultan Drive, Ft. Detrick, MD 21702-5006, DSN 343-7267.

5.2 "Notes" codes in each kit component listing reflect information required by the MMM when performing kit inspections. The most commonly used codes are:

5.2.1 "F" - Item subject to damage by freezing.

5.2.2 "P" - Item with potency period or expiration date.

5.2.3 "Q" - Item which is a drug or other substance requiring security storage.


5.2.4 Other codes may be identified by reviewing the UDR under "Special Notes" prior to the "Table of Contents."

6 SERVICEABILITY TAGS.

6.1 Medical kits listed herein will be tagged with a DD FORM 1574, SERVICEABLE TAG - MATERIEL, or AFTO FORM 104, FIRST AID/SURVIVAL KIT INSPECTION CERTIFICATE.

6.2 The tag or label will be marked as shown in figure 1. Major commands are authorized to grant life support activities the option to delete or obliterate the inspecting activity when life support equipment requires sanitizing. The MMM or the property administrator at an aircraft contractor plant will complete this tag on initial issue of the kit to indicate that the components are up to date and will not expire prior to the next inspection date.

6.3 In addition to inspection information the AFTO FORM 104 is used to protect the integrity of the medical kits. Kits which have a torn AFTO FORM 104 will not be issued by life support or using activities. All such kits will be returned to the MMM for reinspection.

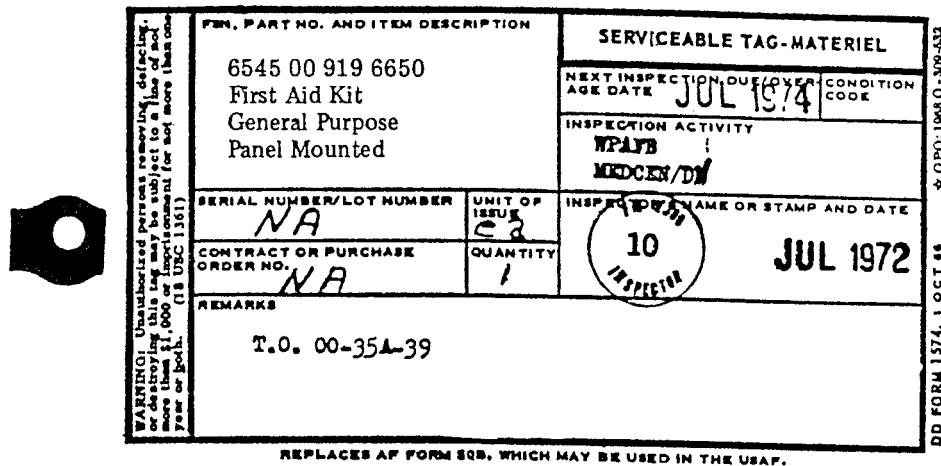

DD FORM 1574, 1 OCT 66

<small>WARNING: Unauthorized persons removing, deleting, or destroying this tag will be subject to court-martial or more than \$1,000 or imprisonment for not more than one year or both. (18 USC 1361)</small>	FORM, PART NO. AND ITEM DESCRIPTION 6545 00 919 6650 First Aid Kit General Purpose Panel Mounted		SERVICEABLE TAG-MATERIEL	
	SERIAL NUMBER/LOT NUMBER NA		UNIT OF ISSUE EA	INSPECTION ACTIVITY WPAFB MEDCEN/DN
	CONTRACT OR PURCHASE ORDER NO. NA		QUANTITY 1	INSPECTOR NAME OR STAMP AND DATE 10 JUL 1972
	REMARKS T.O. 00-35A-39			
	NEXT INSPECTION DUE/ OVERAGE DATE JUL 1974			

REPLACES AF FORM 89B, WHICH MAY BE USED IN THE USAF.

DD FORM 1574, 1 OCT 66

FIRST AID/SURVIVAL KIT INSPECTION CERTIFICATE		
NEXT INSPECTION DUE/ OVERAGE DATE 2 Jan 80		
INSPECTION ACTIVITY Wright-Patterson AFB, Ohio 45433	INSPECTOR FM 2300 10 WPAFB	DATE 2 Jan 1978

AFTO FORM 104
APR 71

H9003172

Figure 1. AFTO Form 104 and DD Form 1574

7 INSPECTION OF MEDICAL KITS.

7.1 Medical kits listed in this TO will normally be inspected by the MMM every 24 months, when suspension of the kit or one of its components is announced in the AFMLL, or when turned in by using activities in accordance with paragraph 8. The 24 month inspection benchmark may be extended or shortened to coincide with aircraft isochronal inspections if sufficient shelf life is available on individual components. If this procedure is followed, inspection of the kit can be included in the aircraft inspection checklist and accomplished during routine aircraft downtime. Medical logistics personnel will make every possible attempt to insure sufficient shelf life is available on kit components under this program. Using activities will insure that kits are removed from use and turned in to the MMM for inspection and quality assurance action prior to the next inspection due and/or overage date or immediately upon receiving notice from the MMM of an applicable item suspension.

7.2 During the inspection the MMM will accomplish the following minimum actions.

7.2.1 Insure kits are properly marked. If not, marking in accordance with paragraph 4 will be accomplished.

7.2.2 Replace all expiration dated items with less than 24 months shelf life remaining wherever possible. If not possible, the NEXT INSPECTION DUE AND/OR OVERAGE DATE annotated on the serviceability tag will be that of the component that has the shortest expiration date or shelf life remaining.

7.2.3 Insure that the kit contains no Schedule 12s controlled medical items unless a current major command authorization is on file.

7.2.4 Insure the contents of the kit are complete and serviceable.

7.2.5 Subsistence items will be replaced by items with a later date of manufacture during normal inspection of the kits.

7.2.6 Insure that a list of contents is included in the kits. Local reproduction of the lists by the MMM is authorized.

7.2.7 Insure that AFR 69-10, Materiel Quality Control Storage Standards, is used as a reference in the inspection of all kits and components.

7.2.8 Component items of first aid and/or survival kits that are deleted through normal supply action will be replaced or withdrawn from the sets and kits at the next normally scheduled inspection cycle. It is not necessary to recall the kits in use to effect replacement of deleted items unless otherwise notified through appropriate channels.

7.2.9 Recommend replacement of the kit if case or all components are unserviceable.

7.2.10 Identify possible abuse of the first aid kit or its contents and advise the appropriate unit commander and the DBMS.

8 INSPECTION BY OTHER THAN MEDICAL PERSONNEL.

Using activities are responsible for periodically inspecting kits in their possession. This inspection is to insure that the kits have not been used, or subjected to an environment hostile to the kit. Kits which have been opened or whose outer appearance indicates possible damage to the contents will be turned in to the MMM for inspection and quality control. Inspections accomplished by using activities will be confined to an external examination of the container and its seals. The DD FORM 1574 or AFTO FORM 104 will be removed only by base medical supply personnel. Inspection actions (paragraph 7) will be coordinated with the MMM to insure a level workload is maintained in both the using activity and medical supply.

9 FUNDING FOR REPLACEMENT COMPONENTS.

9.1 The MMM or his designated representative will replace missing components without charge to requesting USAF activities, including Air Reserve Components. Charges for replacement of component items will be made against RC/CC XX5864, reference DFAS-DER 170-5. The base medical Operations and Maintenance (O&M) program will budget for this cost. If replacement requirements are excessive, the MMM will advise the base commander through the DBMS and request an investigation or other appropriate action.

9.2 The cost of components used to repack kits for other than a USAF activity (including Air Reserve Components) will be charged to the agency requesting the service.

9.3 If exterior cases and/or components are damaged, the requesting USAF activity will fund replacement of the new kit, new case, or affected components.

10 TURN IN OF EXCESS KITS.

10.1 Unserviceable and excess components generated during the inspection by the MMM will be turned in and processed through the stock record account in accordance with AFMAN 23-110, Volume 5, chapter 20.

10.2 Medical first aid kits will be removed from aircraft prior to disposal or transfer of the aircraft to the Military Aircraft Storage and Disposition Center (MASDC). Using activities are responsible for returning these and all other excess or no longer authorized kits to the MMM. Reimbursement to using activities for kits turned in will not be given without the written approval of AFMLO/FOM-P, 1423 Sultan Dr., Ft. Detrick MD 21702-5006.

11 FURNISHING KITS TO AIR FORCE CONTRACTORS.

11.1 Medical kits required for initial installation in new weapons systems will be funded by HQ AFMC. Funded requisitions for kits required by the terms of the Air Force contract will be submitted by the Air Force Plant Representative Office (AFPRO) through HQ AFMC/SGAR to AFMLO/ FOC-T, 1423 Sultan Dr., Ft. Detrick MD 21702-5006, for processing.

11.2 The director of materiel management at each AFMC ALC will compute requirements for kits being installed in aircraft or other major commodities for which they have modernization and/or maintenance responsibilities. These requirements will be funded by HQ USAF/SGHC. The contractor will submit their kit requirement to the director of materiel management at the applicable ALC, who in turn will submit unfunded requisitions to AFMLO/FOC-T, 1423 Sultan Dr., Ft. Detrick MD 21702-5006.

11.3 Contractors having unserviceable kits will request disposition instructions from AFMLO/FOM-P, 1423 Sultan Dr., Ft. Detrick MD 21702-5006.

12 KIT STANDARDIZATION OR MODIFICATION PROCEDURES.

Recommendations to standardize additional medical kits or modify kits will be submitted as prescribed in AFMAN 23-110, Volume 5, chapter 16, paragraph 16-10, or using AFTO FORM 22, TECHNICAL ORDER SYSTEM PUBLICATION IMPROVEMENT REPORT AND REPLY. All such recommendations will be coordinated with both medical and life support authorities at each level of command.

13 NSN 6545-00-116-1410, FIRST AID KIT, GENERAL PURPOSE, RIGID CASE (FIGURE 2).

13.1 **PURPOSE.** An assemblage of items designed to provide first aid material for use in life rafts, airborne lifeboats, large aircraft and small water craft.

13.2 **INSTALLATION.** The USAF activity engaged in assembling survival kits for any of the above designed rafts, watercraft or aircraft is responsible for its installation.

13.3 **MARKING.** This kit will be marked in accordance with paragraph 4 of this TO.

13.4 **INSPECTION.** Using activities will return this kit to the MMM for the inspection required by paragraph 7. The MMM will replace all "P" coded items that have less than two years of shelf life remaining. Follow instructions in paragraph 7. 2 for "P" coded items with less than two years shelf life remaining if they are needed to update the kit.

NOTE

NSN 6545-00-853-6309, first aid kit, eye dressing should be opened and its contents inspected.

13.5 **SEALING AND/OR TAGGING.** Upon completion of the inspection and repacking, the MMM will affix a properly completed AFTO FORM 104 to the kit (figure 2). The AFTO FORM will be affixed in a manner to preclude opening of the kit without tearing the form (figure 2). Kits which have a torn AFTO FORM 104 will not be issued and will be returned to the MMM for reinspection.

14 NSN 6545-00-139-3671, SURVIVAL KIT, INDIVIDUAL (FIGURE 3).

14.1 **PURPOSE.** This kit is designed to provide first aid items, toilet articles, minimal nourishment and notions to individual aircrew members stranded in isolated areas. Special instructions are as follows:

14.2 INSTALLATION.

14.2.1 When used as a supplement to type I survival kit assemblages, the USAF activity engaged in assembling such kits, will install this item in its designated location.

14.2.2 When issued as an item of individual equipment its issue by the individual equipment activity will be recorded on AF FORM 538, PERSONAL CLOTHING AND EQUIPMENT RECORD, or on AF FORM 1297, TEMPORARY ISSUE RECEIPT. Personal equipment activities will establish a method to insure the item is returned every 24 months for inspection and repacking. The kit, when issued as an individual equipment item, will be carried securely in pockets of the outer garment in such a manner that it will not be cumbersome or dangerous to the individual during bail out or ditching procedures.

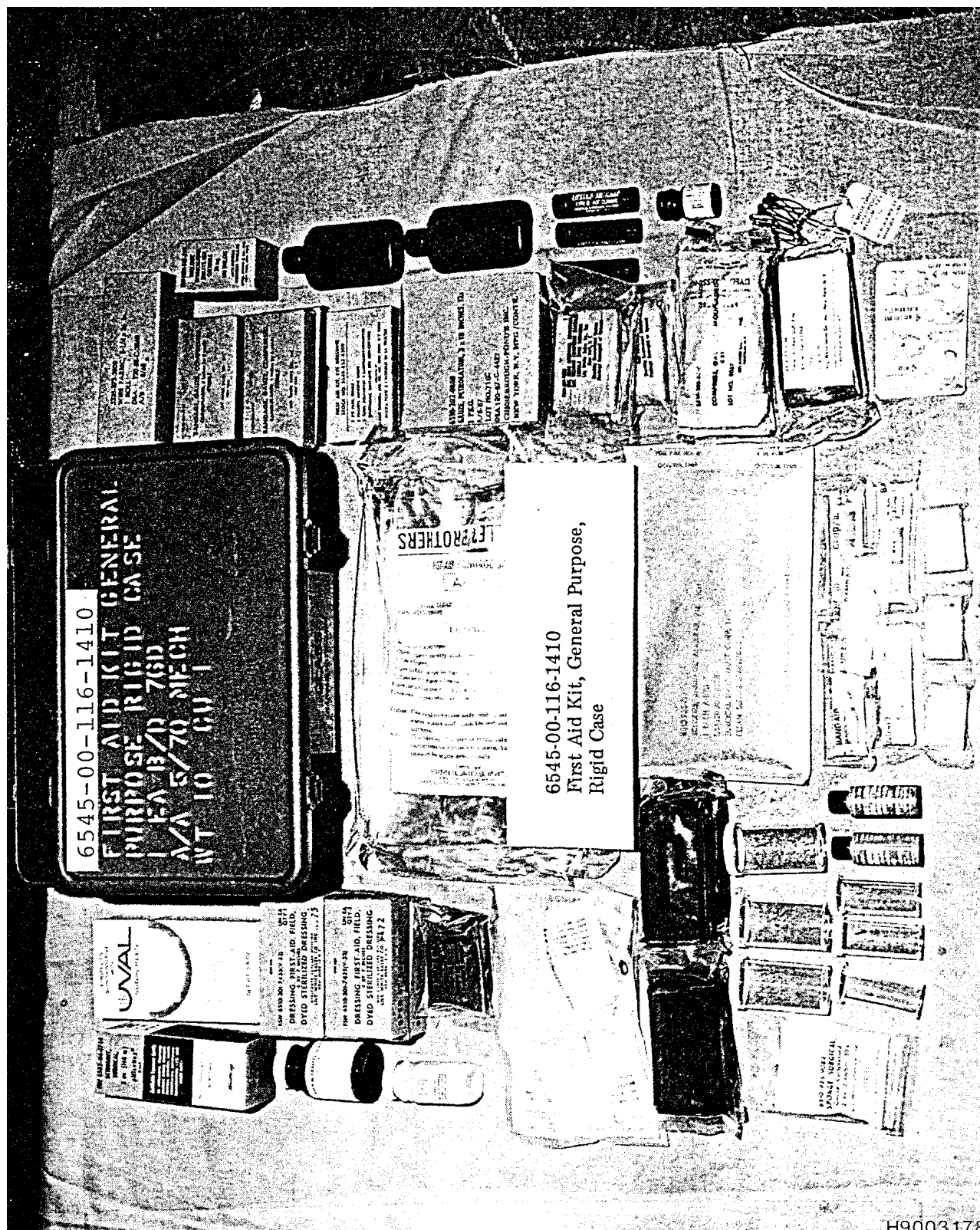


Figure 2. First Aid Kit, General Purpose, Rigid Case, NSN 6545-00-116-1410 (Sheet 2 of 2)

Table 2. NSN 6545-00-116-1410 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-00-131-6992		CASE, MEDICAL INSTRUMENT AND SUPPLY SET, Plastic, Rigid, Size B, 9-1/2 inches long x 9 inches wide x 7 inches high. Can substitute NSN 6545-01-421-7576	1	EA
5335-00-373-2800		WIRE FABRIC, 5-1/4 x 36 inches	1	RO
6505-00-100-9985 (LP)	P	ASPIRIN TABLETS, USP, 0.325 Gram, 100s SOURCE FOR LOCAL PURCHASE: Goldline Laboratories, Inc.; CAGE 63488; 1900 W. Commercial Blvd, Fort Lauderdale FL 33309; Phone: 800-327-4114; Part Number: NDC00182-0444-01, other items available are NDC00395-0134-10, Humco or NDC00904 2009-60, MAJOR.	1	BT
6505-00-491-7557 (LP)	FP	POVIDONE-IODINE CLEANSING SOLUTION, 7-1/2%, 4 fl oz (118 ml). Can substitute NDC 0034-2200-04 SOURCE FOR LOCAL PURCHASE: Purdue Frederick Company; CAGE 05161; 100 Connecticut Avenue, Norwalk CT 06850; Phone: 800-733-1333 or 203-853-0123; Part Number: Betadine Skin Cleanser or NDC 00034-2350-40; must order in multiples of 36 or NDC58879-0135-04, Vitamed Labs.	2	BT
6505-00-926-2111 (LP)	P	MECLIZINE HYDROCHLORIDE TABLETS, USP, Chewable, 25 mg, 100s. SOURCE FOR LOCAL PURCHASE: GOLDLINE, NDC00182-0571-01 OR NDC52152-0117-02, H.L. Moore or Watson/ Rugby NDC00536-3990-01.	1	BT
6505-01-121-2335	P	SUNSCREEN PREPARATION FACTOR 15, 4 fl oz; Part Number NDC24385-0234-26 SOURCE FOR LOCAL PURCHASE: Bergen Brunswick, Schering	2	BT
6508-01-436-0607		CHAPSTICK, Cold or Hot Climate, SFP 15, 4.4 Grams 144s	3/144	PG
6510-00-159-4883		DRESSING, FIRST AID, FIELD, Camouflaged, 4 x 7 inches	4	EA

Table 2. NSN 6545-00-116-1410 - Components - Continued

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6510-00-200-3080		COMPRESS AND BANDAGE, Camouflaged, 4 x 4 inch Compress, 84 inch Bandage	2	EA
6510-00-200-3180		BANDAGE, GAUZE, Compressed, Camouflaged, 2 inches x 6 yards	2	EA
6510-00-200-3190		BANDAGE, GAUZE, Compressed, Camouflaged Olive Green, 4 inches x 6 yards	4	EA
6510-00-201-1755		BANDAGE, MUSLIN, Compressed and Folded, Camouflaged Olive Drab, 37 x 37 x 52 inches, w/Safety Pins	4	EA
6510-00-201-3000		COTTON, PURIFIED, USP, 1 oz, 4-1/2 x 2 inches, Rolled	1	PG
6510-00-201-7680		COMPRESS AND SKULLCAP, HEAD DRESSING	3	EA
6510-00-202-0800		GAUZE, PETROLATUM, 3 x 18 inches, 12s	1	PG
6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	30/300	BX
6510-00-926-9082		SPONGE, SURGICAL, Gauze, Compressed, 2 x 2 inches, 10s	2	PG
6510-01-060-1639	P	ADHESIVE TAPE, SURGICAL, 1 inch x 1-1/2 yards, 100s	3/100	PG
6515-00-687-8052		AIRWAY, PHARYNGEAL, Plastic, Disposable, Adult, 12s	1/12	BX
6515-00-935-6592		SPLINT, ARM, Adult, Pneumatic, Zipper Closure, Nonadjustable, Plastic	1	EA
6515-00-935-6593		SPLINT, LEG, Adult, Pneumatic, Slide Fastener, Nonadjustable, Plastic	1	EA
6515-00-935-7138		SCISSORS, Bandage, Angular, 7-1/4 inch	1	EA
6515-00-958-2232		AIRWAY, PHARYNGEAL, Plastic, Disposable, Child, 12s	1/12	BX
H72678522		GLOVE, PATIENT EXAM/TREAT, Large Size, Prepowdered, Snug Fitting Cuff, 100s	4/100	PG
6520-00-926-9041		BAG, Dental Prosthesis, Plastic, 100s	1/100	PG

Table 2. NSN 6545-00-116-1410 - Components - Continued

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6530-01-316-0282		BOTTLE, Safety Cap, Plastic, 16 drams, Child Resistant, 230s	4/200	PG
6530-00-284-8441		BOTTLE, Safety Cap, Plastic, Light Resistant, Tablet and Capsule, 44 ml, 200s. (The empty plastic bottles in this kit may be used for additional medication as required for operational mission.)	4/200	PG
6545-00-853-6309	P	FIRST AID KIT, EYE DRESSING	1	EA
6840-01-284-3982		INSECT REPELLENT, PERSONAL APPLICATION, 2 oz Tube, 12s	2/12	BX
6850-00-985-7166	P	WATER PURIFICATION TABLET, IODINE, 8 mg, 50s	1	BT
8315-00-787-8000 (GSA Item)		PIN, SAFETY, 1-1/2 inches long, 72s	1	BX
AFTO FORM 401		INSTRUCTION CARD, Artificial Respiration, Mouth-to-Mouth Resuscitation	1	EA
		TOTAL WEIGHT 17.282 TOTAL VOLUME 1.697		

14.3 INSPECTION AND REPLACEMENT OF COMPONENTS. This kit will be inspected and quality assurance actions performed by the MMM in accordance with paragraph 7 of this TO. Using activities will establish procedures to insure all kits are returned to the MMM for this inspection. All items identified with note code "P" will be replaced by the MMM if they have less than two years of shelf life remaining, unless items are not available; follow instructions in paragraph 7.2.2 for "P" coded items with less than two years shelf life remaining if they are needed to update the kit. Subsistence items will be replaced by items containing a later date of manufacture during each inspection of the kits by the MMM.

14.4 SOURCE OF SUPPLY AND SPECIFICATIONS FOR NONSTANDARD ITEMS. Items coded LP in the component list may be procured from the manufacturer of the item or a local distributor. Before purchasing from a local source, check to make sure the item has not been recently depot stocked. Be aware that manufacturers and distributors have minimum order quantities.

14.5 REPACKING OF KIT. All components of this kit are to be enclosed in a waterproof enclosure. Replacement items which are not sealed in a primary plastic or aluminum enclosure will be sealed in polyethylene bags by the MMM. NSN 6850-00-985-7166, Water Purification Tablets, Iodine 8 mg, 50's that are sealed with parafin wax do not require additional sealing. If upon visual inspection, the wax seal is broken, the bottle should be replaced.

Table 3. NSN 6545-00-139-3671 - Components

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
	6545-01-189-7094		CASE, INDIVIDUAL SURVIVAL KIT, Plastic	1	EA
13	4220-01-181-3154		FISHING KIT EMERGENCY (3 Components)	1	EA
11	6505-00-982-0288	P	ERYTHROMYCIN OPHTHALMIC OINTMENT, USP, 5 mg per gram, 1/8 oz, (3.5 gram)	1	TU
3	6510-00-200-3080		COMPRESS AND BANDAGE, CAMOUFLAGED OLIVE GREEN, 4 x 4 inch Compress, 84 inch Bandage	1	EA
12	6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	6/300	BX
	6510-01-060-1639	P	ADHESIVE TAPE, Surgical, 1 inch by 1-1/2 yards, 100s	1/100	PG
	6515-00-296-2529		SAW, Flexible, Steel Wire, 16 Teeth/Inch Long	1	EA
	6515-01-337-3689		GLOVE, Patient Exam/Treatment Mens/Womens Medium Prepowdered Rubber, 100s	4 gloves or 2 pair	PG
7	6850-00-985-7166	P	WATER PURIFICATION TABLETS, IODINE, 8 mg, 50s	1	BT
	7610-01-301-5259 (LP)		LIST, MEDICAL INSTRUCTIONS - CONTENTS	1	SH
2	8520-00-782-3554 (GSA Item)		HAND CLEANER, TOWELETTE, 2000s	3/2000	BX

Table 3. NSN 6545-00-139-3671 - Components - Continued

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
10	8925-01-008-0960 (LP)		<p>CANDY, CHOCOLATE. Candy for the Survival Kit can be any wrapped, high calorie item is acceptable, i.e., granola bar, hard candy, caramels/toffee, fruit bar, brownie, charms, or M&M Mars, Plain M&Ms and Kudo Bars (Granola Bars). Activities can local purchase any suitable item and pack in survival kit. Kudo Bars are available unenrobed (not chocolate covered) and chocolate covered. Unenrobed: Chocolate Chunk, the only other unenrobed bars are Kudo bars in the following flavors: Kudos Milk Chocolate Variety Pack - Part No. 01779 - \$216.00 CS, Kudos Milk Chocolate and Fudge - Part No. 03382 - \$144.00 CS, Kudos Milk Chocolate and Peanut Butter - Part No. 03383 - \$144.00 CS, Kudos Milk Chocolate and Chocolate Chip - Part No. 03384 - \$144.00 CS, Kudos Milk Chocolate w/Mini M&M's - Part No. 03385 - \$144.00 CS. Current Cost is \$0.19 per bar FOB Norfolk, Payment by credit card. Source Opportunity Resources Inc., Candace Lewis, 406-721-2930. New Source of Supply, Specialty Products, Inc., 6206 Easton Rd., Pipersville, PA 18947, POC Fred Prozzillo.</p>	1	BG

Table 3. NSN 6545-00-139-3671 - Components - Continued

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
4	8935-00-082-5667 (LP)		SOUP AND GRAVY BASE, INSTANT, CHICKEN FLAVORED, 8 gram envelopes, 1000s, Managed by DSC-P Subsistence Opportunity Resources, Inc., The Soup and Gravy Base can be ordered from DPSC (S9S) by ordering either of the following NSNs: 8935-01-124-6061, instant, beef flavored, 7 gm or 8935-00-082-5667, instant, chicken flavored, 7 gm. Also the Soup and Gravy Base can be substituted with chicken flavor or, as an alternate, a 5 to 8 gram commercially available individual package of dried broth, bouillon, or soup, suitable for mixing with 6 to 8 ounces of both hot and cold water. If the alternate is used, it shall be packaged in such a manner as to be unaffected by moisture. Packets should be inspected for minute pinholes in foil packs which can cause contamination of contents. Any found to be defective are to be destroyed and replaced with fresh stocks. The soup and gravy base can be purchased from Trans-Packers Service Corp, 419 Vanderfoort Avenue, Brooklyn NY 11222-5313, telephone (718) 963-0900. Their CAGE code is 6D622, do not accept credit cards. The commissary may carry a suitable substitute or be able to order the needed items.	2/1000	PG
5	9535-01-201-7014 (LP)		ALUMINUM FOIL, 12 x 12 inches, Available Opportunity Resources, Inc. POC Candace Lewis, 406-721-2930 and Specialty Products, Inc., POC Fred Prozzillo, 215-766-7777.	1	EA
9	9920-01-191-3434 (LP)		MATCHES, SAFETY, Humidity Resistant, Water Proof Bag, 2 x 1-1/2 inches (1 pkg of 2500 book of matches). The Safety Matches, Humidity Resistant 2500 book of matches, GSA Item can also be purchased from Opportunity Resources Inc., POC Candace Lewis, 406-721-2930 and Dalrex, Inc., 3795 NW 25th Street, Miami, FL 33142, 305-638-8220.	1	PG

Table 3. NSN 6545-00-139-3671 - Components - Continued

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
			TOTAL WEIGHT 1.613 TOTAL VOLUME 0.728		
<p>NSN 4220-01 181-3154, Fishing Kit Emergency (3 Components) Available from Opportunity Resources, POC Candace Lewis, 406-721-2930 and Specialty Products Inc., 6626 Easton Rd. Pipersville, PA 18947, POC Fred Prozzillo at 215-766-777.</p> <p>NSN8935-00-082 5667, Soup and Gravy Base, Instant, Chicken Flavored, 8 gram envelopes 1000s, Managed by DSC-P Subsistence Opportunity Resources, Inc., POC Candace Lewis, 406-721-2930, and TCOMUS Co., POC Bill Chervenak, 931-520-0005.</p>					

14.5.1 The polyethylene bags specified herein shall be made from clear extruded stainless LAYFLAT TUBING. Polyethylene shall be 100 percent virgin material having an average wall thickness of 0.003 inch. Minimum thickness at any point shall be not less than 0.00225 inch. Both seam and closure shall be heat sealed to prevent moisture from entering or to prevent escape of fluid contents in the event of leakage; excess air shall have been excluded from the sealed bag so that bag conforms to the shape of the component therein. To insure a minimum of excess material, the width of commercially available tubing shall be of a dimension which shall closely fit the individual component to be enclosed. The length of the tubing between seals shall be the minimum to contain the component. Excess tubing beyond the seals shall be cut off. A properly sealed bag can be submerged to a depth of one foot in tap water at room temperature for 15 minutes without leaking.

14.5.2 Attach tape handles to kit by using tape pressure sensitive adhesive, adhesive filament-reinforced, 3/4 inch wide, NSN 7510-00-290-8035 low tensile strength or NSN 7510-00-802-8311 medium tensile strength. If the 3/4 inch wide tape is not available, you may use one inch wide NSN 7510-00-290-8036 low tensile strength or a NSN with medium tensile strength. These handles will be added lengthwise to the bottom of the case and another attached to the top cover at right angles as shown in figure 3. Handles are added to provide a carrying handle, for ease in opening, and for attachment of kit to survival vest mesh net.

14.6 MARKING. This kit does not require a red cross symbol as it is an individual item of equipment.

14.7 SEALING AND/OR TAGGING. Upon completion of the inspection and repacking the MMM will affix a properly completed AFTO FORM 104 to the kit. The AFTO FORM will be affixed in a manner to preclude opening of the kit without tearing of the form. Kits which have a torn AFTO FORM 104 will not be issued. All such kits will be returned to the MMM for reinspection.

15 NSN 6545-00-782-6412, SURVIVAL KIT, INDIVIDUAL, TROPICAL, TACTICAL AIRCREW (FIGURE 4).

15.1 PURPOSE. This kit is specifically designed to provide a one week supply of first aid and survival supplies to stranded aircrew members during escape and evasion situations encountered under tropical conditions.

15.2 INSTALLATION. This kit will be issued in accordance with paragraph 14.2 of this TO.

15.3 INSPECTION. Using activities will return this kit to the MMM for inspection required by paragraph 7. The MMM will replace all "P" coded items that have less than two years of shelf life remaining unless items are not available. Follow instructions in paragraph 7.2 for "P" coded items with less than two years of shelf life remaining if they are needed to update the kit.

15.4 REPACKING.


15.4.1 Components of this kit are installed as shown in figure 4. Aluminum foil strip items are to be folded to fit the designated places. Rubber bands are to be used to keep the strips folded. Items which are not encased in aluminum foil will be encased in a polyethylene bag as prescribed in paragraph 14.5.1.

15.4.2 The inner surface of the bottom piece of the plastic case is coated with a layer of non-curing, pressure sensitive adhesive. This adhesive is provided to secure positioning of the components in their respective locations while the case is opened or inverted.

15.4.3 Items coded LP in the component listing may be procured from the manufacturer of the item or a local distributor. Before purchasing from a local source, check to make sure the item has not been recently depot stocked. Be aware that manufacturers and distributors have minimum order quantities.

15.5 MARKING. This kit does not require application of a red cross symbol.

15.6 SEALING AND/OR TAGGING. After completion of the inspection and repacking, the kit will be tagged with two AFTO FORMs 104 - one on each side of the kit, see figure 1. The AFTO FORMs 104 will be affixed in a manner to preclude opening the kit without tearing the forms. Kits which have a torn AFTO FORM 104 will not be issued by life support or using activities. All such kits will be returned to the MMM for reinspection.

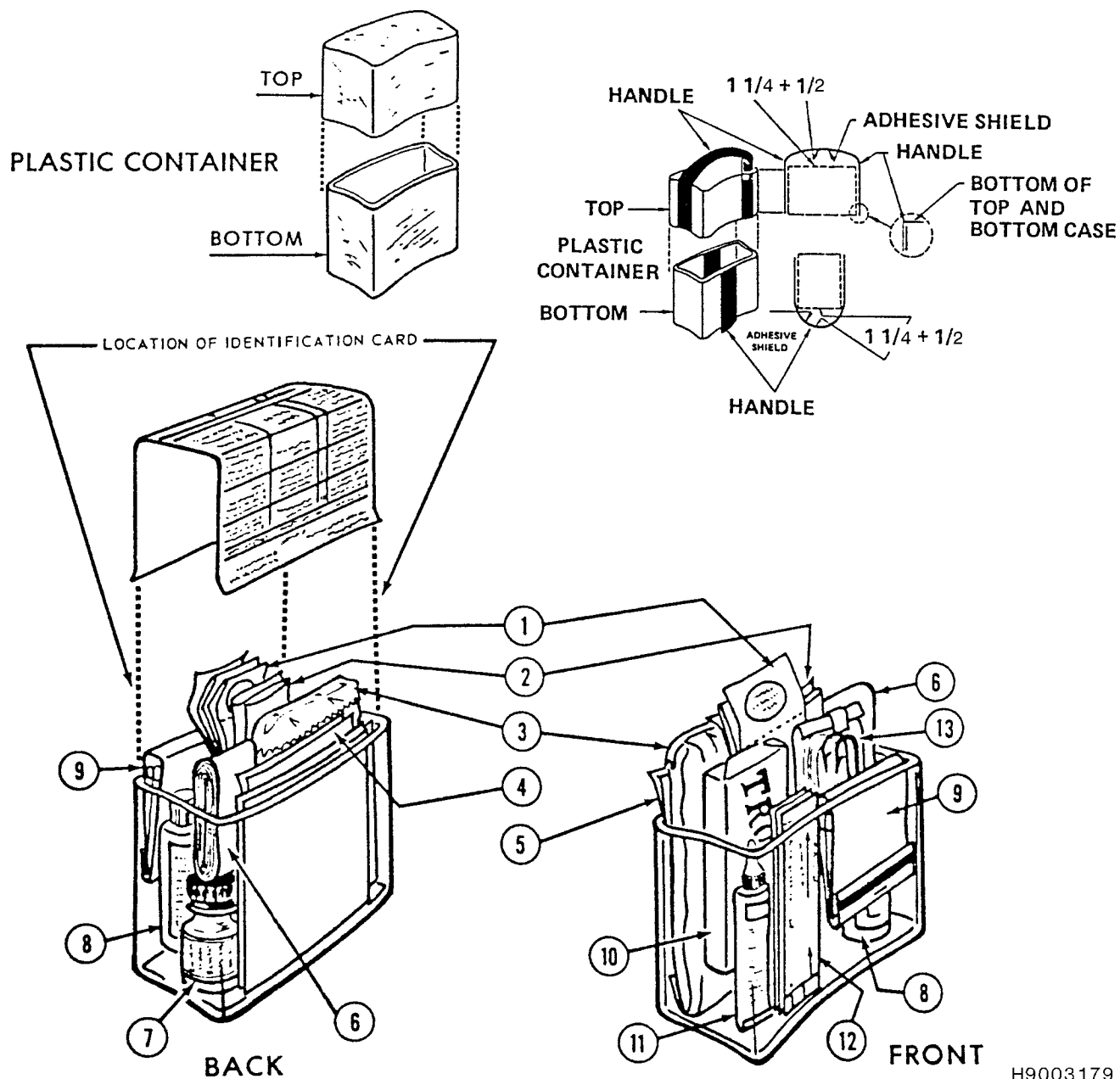

H9003179


Figure 3. Kit Assembly

Table 4. NSN 6545-00-782-6412 - Components

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
	6545-01-189-7095		CASE, INDIVIDUAL, Survival Kit, Plastic. New Source of Supply, Specialty Products, Inc., 6626 Easton Rd., Pipersville, PA 18947, POC Fred Prozzillo, 215-766-7777	1	EA
	4220-01-181-3154		FISHING KIT EMERGENCY (3 Components). Available from Opportunity Resources, POC Candace Lewis, 406-721-2930 and Specialty Products Inc., 6626 Easton Rd., POC Fred Prozzillo, 215-766-7777	1	EA
3	6505-00-009-5060 (LP)	P	DOXYCYCLINE HYCLATE CAPSULES, USP, 100 mg, I.S., 100s (Anti-Malarial Prophylactic Drug)	14/100	PG
6b	6505-01-238-5632 (LP)	P	LOPERAMIDE HYDROCHLORIDE CAPSULE, 2.0 milligrams Loperamide, I.S., 100s (Anti-Diarrhea)	10/100	PG
8	6505-00-982-0288	P	ERYTHROMYCIN OPHTHALMIC OINTMENT, USP, 5 mg per gram, 1/8 oz (3.5 gram)	1	TU
12	6510-00-913-7906		BANDAGE, GAUZE, ELASTIC, Two-Ply, 2 inches by 5 yards, 50s	1/50	BX
13	6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 by 3 inches, 300s	12/300	BX
14	6510-01-060-1639	P	ADHESIVE TAPE, SURGICAL, 1 inch by 1-1/2 yards, 100s	1/100	PG
9	6515-00-296-2529		SAW, Flexible, Steel Wire Cadmium Plated, 16 teeth, 16 Inch Long	1	EA
11	6515-01-363-1212		RAZOR, SURGICAL PREPARATION, Small, Folding Handle, Straight Type Blade, 100s	1/100	PG
	6515-01-337-3689		GLOVE, PATIENT EXAM/TREATMENT, mens/womens, Medium, Prepowdered, rubber, 100s	4 gloves or 2 pair	PG

Table 4. NSN 6545-00-782-6412 - Components - Continued

ITEM NO.	NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
10	6605-01-326-1654		COMPASS, Magnetic, Unmounted, Pocket Type, Source of Supply - Stocker and Yale, Inc., J-483 000, 603-893-8778 and Commega Association, telephone 616-392-7999, Part Number J483.	1	EA
21	6695-00-118-2918 (LP)		BAG, FOOD SAMPLE, Polyethylene, Disposable, 18 fl oz Capacity, 500s	5/500	PG
17	6840-00-142-8965		INSECT REPELLENT, Personal Application, 1 oz, 12s	1/12	BX
5	6850-00-985-7166	P	WATER PURIFICATION TABLET, IODINE, 8 mg, 50s	1	BT
19	7610-01-301-5257		LIST MEDICAL INSTRUCTIONS - CONTENTS (Local Reproduction Authorized)	1	SH
16	8415-01-192-2357		INSECT NET MITTENS SET, Head and Hand Covering, Nylon, Black, Source of Supply - Specialty Products Inc., 6626 Easton Rd., Pipersville, PA 18947, POC Fred Prozzillo, 215-766-7777	1	EA
20	9920-01-191-3434 (LP)		MATCHES, Safety, Humidity Resistant, Waterproof Bag 2x1-1/2 inches (2500) Source of Supply: GSA and Opportunity Resources Inc., POC Candace Lewis, 406-721-2930 and Datrex Inc., 3795 NW 25th Street, Miami, FL, 33142, 305-638-8220. TOTAL WEIGHT 1.069 TOTAL VOLUME 0.341	1	PG

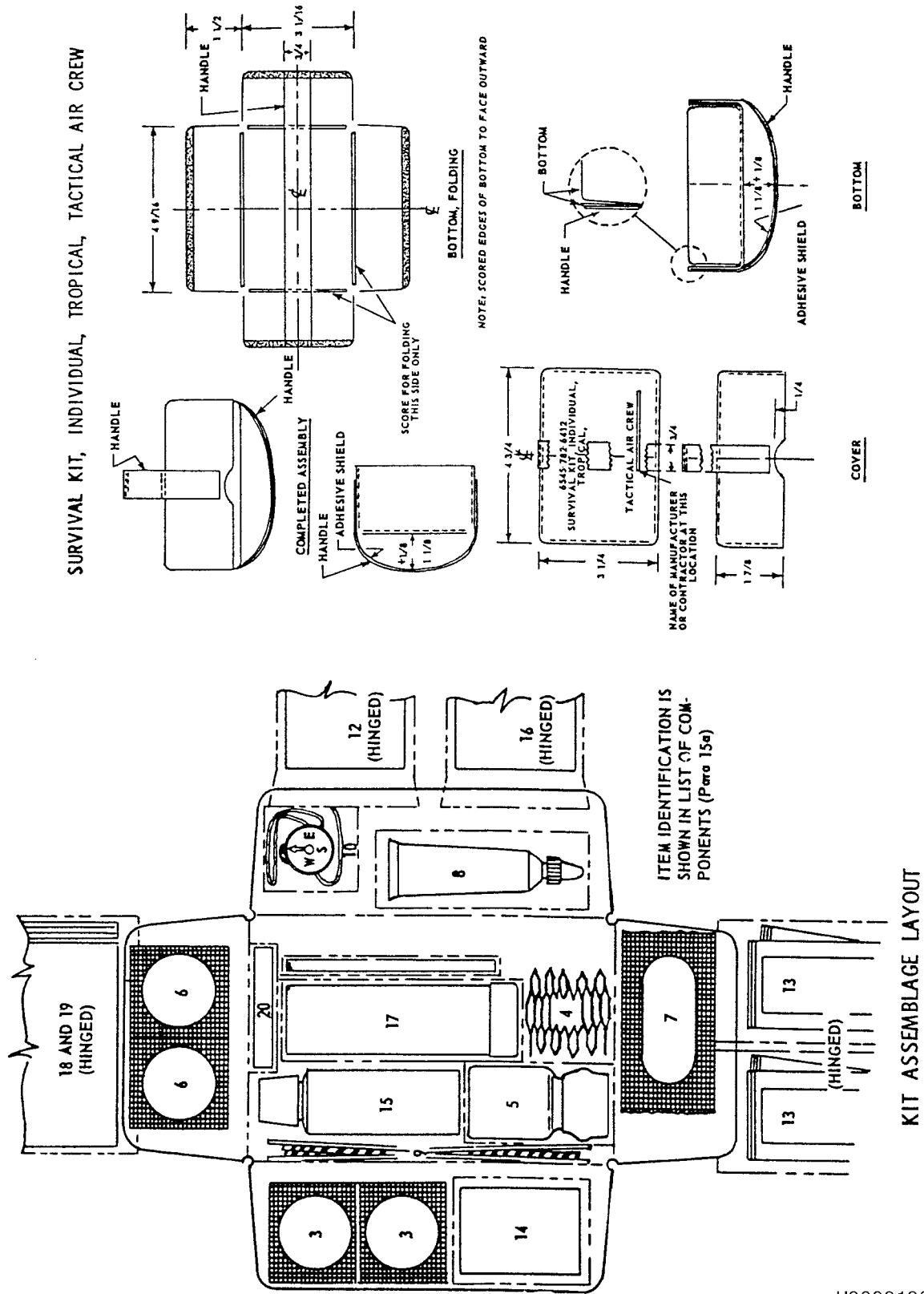

H9003180

Figure 4. Survival Kit, Individual, Tropical, Tactical Aircrew, NSN 6545-00-782-6412

Table 5. NSN 6545-00-853-6309 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6505-00-982-0288	P	ERYTHROMYCIN OPHTHALMIC OINTMENT, USP, 5 mg per gram, 1/8 oz (3.5 gram)	2	TU
6510-01-060-1639	P	ADHESIVE TAPE, Surgical, 1 inch wide by 1.5 yards long, 100s	2/100	PG
6510-01-107-7575		PAD, Cotton Gauze Covered, 2.5 by 2-1/8 inches, 50s	4/50	PG
8105-00-837-7754 (GSA Item)		BAG, Plastic, Interlock Seal Closure Fastening, 6 by 6 inches	1/1000	MX
		TOTAL WEIGHT 0.151 TOTAL VOLUME 0.008		

16 NSN 6545-00-853-6309, FIRST AID KIT, EYE DRESSING.

16.1 PURPOSE. This kit is designed to provide minimal first aid supplies to treat eye injuries. It is not normally issued as an end item but is a component of other kits listed in this TO.

16.2 INSTALLATION. The MMM will install this item in first aid kits when it is listed as a component of the major assembly.

16.3 INSPECTION. The kit will be inspected using the criteria of paragraph 7 whenever the major assembly is inspected.

16.4 MARKING. No special markings are required.

16.5 SEALING AND/OR TAGGING. This item does not require tagging.

16.6 SPECIAL INSTRUCTIONS. The eye dressing kit contains four items, one of which is potency dated. When received from the DPSC, the kit is sealed in a zip lock bag with an additional perforated seal. When servicing the eye dressing kit you must tear the perforated seal and open the zip lock. At that point, simply inspect and/or replace two tubes of erythromycin ophthalmic ointment, USP, 5 mg per gram, 1/8 oz (3.5 gram), NSN 6505-00-982-0288 and reseal the zip lock. There is no need for additional marking and/or sealing of the kit. Each time the prime kit is inspected, ensure sufficient dating is available on the ointment.

17 NSN 6545-00-919-6650, FIRST AID KIT GENERAL PURPOSE, FLEXIBLE CASE, PANEL MOUNTED (FIGURE 6).

17.1 PURPOSE. This kit is an assemblage of items designed to provide first aid supplies to personnel on board an aircraft or watercraft.

NOTE

This item is not suitable for being carried on the person during bail out procedures.

17.2 INSTALLATION. The USAF activity responsible for equipping waterborne craft or aircraft with authorized first aid and/or survival equipment will install this kit in accordance with applicable TOs, parts lists or drawings. Parts or material required to provide brackets or holders for installation of the kit will be obtained from Air Force stock or locally manufactured. The maintenance or repair activity will install the mountings.

17.3 INSPECTION. This kit will be inspected in accordance with paragraph 7.

17.4 REPACKING. The contents of this kit will be installed in the sections indicated in the component list.

17.5 TAGGING. Prior to issue of this kit, it will be tagged by medical material personnel as follows:

17.5.1 Close the zipper completely with the side fastener tab in the down and locked position.

17.5.2 Using seal, NSN 5340-01-004-9387, pass through eyelet of DD FORM 1574 and through the hole of the zipper slide fastener, then lock.

17.5.3 If the operation is correctly performed you will be unable to open the zippered section of the kit without breaking the seal. This method provides positive external visual indicator of kit serviceability.

17.6 MARKING. The outer container of the kit will be marked in accordance with paragraph 4.

Table 6. NSN 6545-00-919-6650 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-00-912-9860		CASE, MEDICAL, INSTRUMENT AND SUPPLY SET, Nylon, Nonrigid, No. 2, 7-1/2 inches long x 4-3/8 inches wide x 4-1/2 inches high	1	EA
IN UPPER POCKET				
6510-00-200-3075		COMPRESS AND BANDAGE, Camouflaged, Compress 2 x 2 inches, Bandage 48 x 2 inches, 4s	1	PG
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, Camouflaged, 37 x 37 x 52 inches, w/Safety Pins	1	EA
6510-01-112-6414		GAUZE, PETROLATUM, 3 x 36 inches, 12s	3/12	PG
6515-00-754-0426		BLADE, SURGICAL PREPARATION RAZOR, Straight Single Edge, 5s	1	PG
UPN H72678522		GLOVE, Patient Exam/Treatment, Large Size, Prepowdered, Snug Fitting Cuff, 100s. Source of Supply: Universal Home Health and Industrial Supplies, Inc., 333 Fallenburg Rd., Tampa, FL, 33169, 813-661-2888.	4/100	PG
6520-00-926-9041		BAG, Dental Prosthesis, Plastic, 100s	1/100	PG
IN LOWER POCKET				
6510-00-159-4883		DRESSING, FIRST AID, FIELD, Camouflaged, 4 x 7 inches	3	EA
6510-00-200-3185		BANDAGE, GAUZE, Compressed, Camouflaged, 3 inches x 6 yards	2	EA
6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	18/300	BX
6510-01-477-9736	P	PAD, POVIDONE-IODINE IMPREGNATED, 2 1/2 in X 1 1/2 inch, 100s,	10/100	PG
6510-01-060-1639	P	ADHESIVE TAPE, Surgical, 1 inch x 1-1/2 yards, 100s	3/100	PG

Table 6. NSN 6545-00-919-6650 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-00-853-6309	P	FIRST AID KIT, EYE DRESSING	1	EA
AFTO FORM 401		INSTRUCTION CARD, Artificial Respiration (Mouth-to-Mouth Resuscitation and First Aid in English)	1	EA
		INSTRUCTION SHEET AND LIST OF CONTENTS (Local Reproduction Authorized)	1	EA
		TOTAL WEIGHT 1.919 TOTAL VOLUME 0.115		



Figure 5. First Aid Kit, General Purpose, Flexible Case, Panel Mounted, NSN 6545-00-919-6650 (Sheet 1 of 2)



Figure 5. First Aid Kit, General Purpose, Flexible Case, Panel Mounted, NSN 6545-00-919-6650 (Sheet 2 of 2)

18 NSN 6545-00-922-1200, FIRST AID KIT GENERAL PURPOSE, RIGID CASE (FIGURE 7).

18.1 PURPOSE. This kit is an assemblage designed to provide a source of emergency first aid items for vehicle operators or passengers, and to satisfy requirements specified for use in AFOSH Standard 127-50, Ground Communications - Electronics (C-2) Systems, paragraph 4-f(5)(a)4.

18.2 INSTALLATION. This kit will be installed by the local maintenance people in a manner that will assure that its container will be firmly retained by locally provided brackets but which will allow the kit to be readily removed without tools for inspection or an emergency.

18.3 INSPECTION. This kit will be returned every 24 months to the MMM for inspection as outlined in paragraph 7.

18.4 MARKING. This kit will be marked as outlined in paragraph 4.

18.5 SEALING AND/OR TAGGING. After completion of the inspection and repacking, the kit will be tagged with AFTO FORM 104 (figure 1). The AFTO FORM 104 will be affixed in a manner to preclude opening the kit without tearing the form. Kits which have a torn AFTO FORM 104 will not be issued by life support or using activities. All such kits will be returned to MMM for reinspection.

Table 7. NSN 6545-00-922-1200 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-00-113-3722 (LP)		CASE, MEDICAL INSTRUMENT AND SUPPLY SET, Plastic, Size A, 7.5 x 4.5 x 2.75 inches, Olive Drab	1	EA
6505-00-148-7096	P	POVIDONE-IODINE OINTMENT, USP, 10%, 1/8 oz, (3.54 grams), 144s	10/144	PG
6505-00-159-4883		DRESSING, FIRST AID FIELD, Camouflaged, 4 x 6-1/4 to 7-1/4 inches	3	EA
6510-00-200-3075		COMPRESS AND BANDAGE, Camouflaged, Compress 2 x 2 inches, Bandage 48 inches, 4s	1	PG
6510-00-200-3185		BANDAGE GAUZE, Compressed, Camouflaged, 3 inches x 6 yards	2	EA
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, Camouflaged Brown, 37 x 37 x 52 inches, with Safety Pins	2	EA
6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	18/300	BX
6510-00-935-5820		BANDAGE, ELASTIC, Cotton, Rolled, 2 inches x 4.5 yards, 12s	1/12	PG
6510-01-060-1639	P	ADHESIVE TAPE, Surgical, 1 inch x 1-1/2 yards, 100s	3/100	PG
6510-01-112-6414		GAUZE, PETROLATUM, 3 x 36 inches, 12s	3/12	PG
6515-00-754-0426		BLADE, SURGICAL PREPARATION RAZOR STRAIGHT, Single Edge, 5s	1	PG
6515-00-935-7138		SCISSORS, Bandage, 7-1/4 inches long	1	EA
UPN H72678522		GLOVE, Patient Exam/Treat, Large Size, Prepowdered, Snug Fitting Cuff, 100s. Source of Supply: Universal Home Health and Industrial Supplies, Inc., 333 Fallenburg Rd., Tampa, FL 33619, 813-661-2888.	4/100	PG
6520-00-926-9041		BAG, Dental Prosthesis, 100s	1/100	PG
6545-00-853-6309	P	FIRST AID KIT, EYE DRESSING	1	EA
AFTO FORM 401		INSTRUCTION CARD, Artificial Respiration (Mouth-to-Mouth Resuscitation and First Aid in English)	1	EA
NONE		LIST OF CONTENTS, (Local Reproduction Authorized)	1	EA
		TOTAL WEIGHT 4.593 TOTAL VOLUME 0.247		



Figure 6. First Aid Kit, General Purpose, Rigid Case, NSN 6545-00-922-1200 (Sheet 1 of 2)



Figure 6. First Aid Kit, General Purpose, Rigid Case, NSN 6545-00-922-1200 (Sheet 2 of 2)

19 NSN 6545-00-927-3000, MEDICAL INSTRUMENT AND SUPPLY SET, AERONAUTIC EMERGENCY.

19.1 **PURPOSE.** This kit is a compact assemblage of sophisticated medical supplies suitable for aerial delivery to, or use by, aircrew members who are forced to bail out.

19.2 **INSTALLATION.** The USAF activity responsible for assembling Global Survival Kits, Aerial Delivery Survival Kits, or Rescue Boat Emergency Equipment will install this kit.

19.3 **INSPECTION.** This item will be inspected every 24 months as prescribed in paragraph 7.

19.4 **MARKING.** This kit will be marked in accordance with paragraph 4.

19.5 **SEALING AND/OR TAGGING.** This kit will be tagged after inspection with DD FORM 1574 or AFTO FORM 104. The tag will be attached by means of a wire loop.

20 NSN 6545-01-090-0645, FIRST AID SET, SHELTER.

20.1 **PURPOSE.** This kit provides an assembly of first aid shelter supplies for field use.

20.2 **INSTALLATION.**

20.2.1 Materiel requirements for this purpose will consist of one set, first aid, shelter, NSN 6545-01-090-0645 and six each litter, rigid pole, for each 100 programmed military personnel.

20.2.2 Items prepositioned under this program will be maintained by the supporting stock record account. The major commands will develop arrangements to insure materiel availability to satellite and collocated operating base personnel.

Table 8. NSN 6545-00-927-3000 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-00-912-9830		CASE, MEDICAL INSTRUMENT AND SUPPLY SET, Plastic	1	EA
6505-00-254-5527 (LP)	P	PETROLATUM, WHITE, USP, 1 oz	1	TU
6505-01-454-2363	P	POVIDONE-IODINE TOPICAL SOLUTION, USP, 10%, 15ml, BT	3	BT
6510-00-159-4883		DRESSING, FIRST AID, FIELD, Individual, Camouflaged, 4 x 6-1/4 x 7-1/4 inches	4	EA
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, Camouflaged, Brown, 37 x 37 x 52 inches, with Safety Pins	2	EA
6510-00-201-7430		DRESSING, FIRST AID, FIELD, Camouflaged, 7-3/4 x 7-1/4 inches, with Safety Pins	2	EA
6510-00-203-5000	P	ADHESIVE TAPE, SURGICAL, 3 inches x 5 yards	1	SL
6515-00-333-2200 (LP)		FORCEPS, DRESSING, 3-1/2 inch	1	EA
6515-00-373-2100		SPLINT, WIRE, Aluminum, 3-3/4 x 30 inches	2	EA
6515-00-383-0565		TOURNIQUET, NON-PNEUMATIC, Camouflaged, 1-1/2 x 41-1/2 inches	1	EA
6515-00-687-8052		AIRWAY, PHARYNGEAL, Adult, Plastic, Disposable, 12s	1/12	BX
6515-00-754-0426		BLADE, SURGICAL PREPARATION RAZOR, Straight, Single Edge, 5s	1	PG
UPN H72678522		GLOVE, PATIENT EXAM/TREATMENT, Large Size, Prepowdered, Snug Fitting Cuff, 100s, Source of Supply: Universal Home Health and Industrial Supplies, Inc., 333 Fallenburg Rd., Tampa FL 33619, 813-661-2888.	4/100	PG
6520-00-926-9041		BAG, Dental, Prosthesis, 100s	1/100	PG
6545-00-853-6309	P	FIRST AID KIT, EYE DRESSING	2	EA
6850-00-985-7166	P	WATER PURIFICATION TABLET, Iodine, 8 mg, 50s	1	BT
8105-01-195-8730 (GSA Item)		BAG, PLASTIC, 24 x 24 inches, 100s. GSA, A-A-1668.	1/100	BX
8315-00-787-7000 (GSA Item)		PIN, SAFETY, Brass-Nickel, 2 inches, 48s	12/48	BX
		TOTAL WEIGHT 4.283 TOTAL VOLUME 0.201		

Table 9. NSN 6545-01-090-0645 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6505-01-454-2363	P	POVIDONE-IODINE TOPICAL SOLUTION, USP, 10%, 15ml, BT	1	BT
6510-00-111-0708		PAD, Non-Adherent, 3 by 4 inches, 100s	1	PG
6510-00-159-4883		DRESSING, FIRST AID, Field, Camouflaged, 4 x 7 inches	10	EA
6510-00-200-3075		COMPRESS AND BANDAGE, Camouflaged, Compress 2 x 2 inches, Bandage 48 x 2 inches, 4s	3	PG
6510-00-200-5000		BANDAGE, GAUZE, Rolled, 3 inches x 10 yards, 12s	1	PG
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, 37 x 37 x 52 inches, with Safety Pins	10	EA
6510-00-201-7430		DRESSING, FIRST AID, Field, Camouflaged, 7-1/2 x 8 inches	10	EA
6510-00-203-5500	P	ADHESIVE TAPE, Surgical, 12 inches x 10 yards	1	RO
6510-00-597-7469	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 100s	1	BX
6510-00-782-2698		SPONGE, SURGICAL, Gauze, Compressed, 4 x 4 inches, 200s	1	PG
6510-00-935-5822		BANDAGE, COTTON, Elastic, Rolled, 4 inches x 4-1/2 yards, 12s	1	PG
6515-00-383-0565		TOURNIQUET, Non-Pneumatic, Camouflaged, 1-1/2 x 42 inches	2	EA
6515-00-935-6593		SPLINT, LEG, Pneumatic, Adult, Slide Fastener	5	EA
6515-00-935-7138		SCISSORS, Bandage, Angular, 7-1/4 inch	5	EA
6515-01-225-4681		SPLINT, UNIVERSAL, Aluminum, 36 x 4-1/2 inches, Light Weight, 12s	1	PG
UPN H72678522		GLOVE, Exam N/S Latex Bisque Large	10/100	PG
6520-00-926-9041		BAG, Dental Prosthesis, 100s	1/100	PG
6530-00-783-7905		LITTER, Folding, Rigid Pole	6	EA
6530-00-784-4205		STRAP, Webbing, Patient Securing, Buckle with Locking Device and Spring, Olive Drab	12	EA
6545-00-914-3500		CHEST, Misc Field, Med Instrument and Supply Set, No. 5, 30 x 18 x 16 inches	1	EA
7210-00-715-7985		BLANKET, Bed, Wool, Olive Green	6	EA
8315-00-787-8000 (GSA Item) AFTO FORM 401		PIN, Safety, 1-1/2 inch, 72s	24/72	BX
		INSTRUCTION CARD, Artificial Respiration, Mouth-to-Mouth Resuscitation and First Aid	1	EA

Table 9. NSN 6545-01-090-0645 - Components - Continued

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
		TOTAL WEIGHT 115.021		
		TOTAL VOLUME 12.863		
NOTE: NSN 6530-00-783-7905 in accordance with AFMAN 23-110, Vol 5.6.3				

20.2.3 This materiel will be issued to supported units only under predetermined defense conditions as specified by the major command. Accordingly, local plans will include procedures for rapid issue of this materiel upon direction of the major command.

20.3 INSPECTION. This item will be inspected every 24 months as prescribed in paragraph 7.

20.4 MARKING. This kit will be marked in accordance with paragraph 4.

20.5 SEALING AND/OR TAGGING. Upon completion of the inspection and repacking the MMM will affix a properly completed AFTO FORM 104 to the kit. The AFTO FORM will be affixed in a manner to preclude opening of the kit without tearing of the form. Kits which have torn AFTO FORM 104 will not be issued. All such kits will be returned to the MMM for reinspection.

21 NSN 6545-01-094-8412, FIRST AID KIT, INDIVIDUAL.

Packed in case, 6545-01-094-6142.

21.1 PURPOSE.

21.1.1 The kit provides first aid items for aircrew members to treat injuries sustained while separated from their aircraft during emergencies. Preferred kit has been enlarged for sufficient space to hold all components.

21.1.2 Kit is a component of standard mobility bag.

21.1.3 This kit is to be included as the individual first aid kit for mobility bags. NSN 6545-01-094-8412 is considered complete for the basic Self Aid and Buddy Care (SABC) concept.

21.2 INSTALLATION.

21.2.1 When used in various major survival kits, or standard mobility bags, the USAF activity responsible for assembling the major items will install this kit.

21.2.2 When issued as an item of individual equipment on AF FORM 538 or 1297, the kit will be carried in an outer pocket or on a vertical strap of the field clothing except during bail out procedures.

21.3 INSPECTION. Using activities will return this kit to the MMM for inspection required by paragraph 7. The MMM will replace all "P"coded items with less than two years of shelf life remaining unless items are not available. Follow instructions in paragraph 7.2 for "P"coded items with less than two years shelf life remaining if they are needed to update the kit.

21.4 MARKING. No special markings are required.

21.5 PACKING AND REPACKING. Components of this kit are to be packed in the plastic insert (6545-01-094-6136). If all contents do not fit, place the plastic insert and the overflow in polyethylene reusable interlocking seal closure bag 8105-00-837-7755.

21.6 SEALING AND/OR TAGGING. After completion of the inspection and repacking, AFTO FORM 104 will be used to tag the kit. The AFTO FORM 104 will be affixed NSN 6545-01-094-6136 insert, first aid kit case, plastic, olive drab, in a manner to preclude opening the kit without tearing the form, or, the kit will be tagged with AFTO FORM 104 across the opening of the polyethylene bag. Kits which have a torn AFTO FORM 104 will not be issued by life support or using activities. All such kits will be returned to the MMM for reinspection.

Table 10. NSN 6545-01-094-8412 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-01-094-6142		CASE, MEDICAL INSTRUMENT AND SUPPLY SET, Nylon, Nonrigid, No. 8, 4-3/4 inches long x 2-3/8 inches wide x 5-1/2 inches high	1	EA
6545-01-094-6136		INSERT, FIRST AID KIT CASE, Plastic, Olive Drab, 4-3/4 inches long x 2-3/8 inches wide x 5-1/2 inches high, with attached cover	1	EA
6505-01-454-2363	P	POVIDONE-IODINE TOPICAL SOLUTION, USP, 10%, 15ml, BT	1	BT
6508-01-436-0607		LIPSTICK, ANTICHAP, Cold or Hot Climate, SPF 15, 4.4 Gram, 144s	1/144	PG
6510-00-159-4883		DRESSING, FIRST AID, Field, Camouflaged, 4 x 7 inches	2	EA
6510-00-200-3180		BANDAGE, GAUZE, Compressed and Folded, Brown, 2 inches x 6 yards	1	EA
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, Brown, Camouflaged, 37 x 37 x 52 inches, Triangular w/Safety Pins	1	EA
6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	18/300	BX
6530-01-316-0282		BOTTLE, Safety Cap, Plastic, 16 Drams, 230s. Source of Supply: Kerr Group Inc., Distributor Don and Richards, 800-722-4866, Part Number CR-16, Cage OM7F6, \$33.00 CS of 230.	2/200	PG
6545-00-853-6309	P	FIRST AID KIT, EYE DRESSING	1	EA
6850-00-985-7166	P	WATER PURIFICATION TABLET, IODINE, 8 mg, 50s	1	BT
AFTO FORM 401		INSTRUCTION CARD, Artificial Respiration, Mouth-to-Mouth Resuscitation	1	EA
		TOTAL WEIGHT 1.285 TOTAL VOLUME 0.100		
NOTE: The empty plastic bottles in this kit may be used for additional medications as required for operational missions.				

22 NSN 6545-01-400-3397, FIRST AID KIT, INDIVIDUAL, MODIFIED.

Packed in case, 6545-01-094-6142.

22.1 PURPOSE

22.1.1 The kit, when completed with dated items and antichap lipstick components of kit 6545-01-094-8412, provides first aid items for aircrew members to treat injuries sustained while separated from their aircraft during emergencies.

22.1.2 The kit is a modification of kit 6545-01-094-8412 to remove dated items and antichap lipstick. The intent of this action is to reduce manhours required for frequent inspection of dated items in the original kit as well as eliminate waste of expired stock. At the time of mobilization, the dated items shall be issued along with the BW/CW antidotes. If the deployment is to an area not requiring BW/CW, the dated items are still required to be issued to complete the original kit. The MMM is responsible for maintaining and issuing the dated items and ensuring that kit 6545-01-094-8412 is the deploying kit.

22.1.3 This modified kit is to be maintained and inventoried with components shown in Table 11 only. Upon issue, the MMM will ensure the dated items and antichap lipstick are added to become the deployed kit 6545-01-094-8412, the individual first aid kit for mobility bags. It is then considered complete for the basic Self Aid and Buddy Care (SSABC) concept.

22.2 INSTALLATION. This kit is not installed in its modified form. The completed kit, 6545-01-094-8412, is the issuable kit for major survival kits or standard mobility bags. Follow installation instructions for this kit.

22.3 INSPECTION. This modified kit requires inspection by the MMM every 24 months as required by paragraph 7.

22.4 MARKING. No special markings are required.

22.5 PACKING AND REPACKING. Components of this kit are to be packed in the plastic insert (6545-01-094-6136). When dated items are added for issue of the kit, if all contents do not fit, place the plastic insert and the overflow in polyethylene reusable interlocking seal closure bag 8105-00-837-7755.

22.6 SEALING AND/OR TAGGING. After completion of the inspection and repacking, AFTO FORM 104 will be used to tag the kit. The AFTO FORM 104 will be affixed to NSN 6545-01-094-6136 insert, first aid kit case, in a manner to preclude opening the kit without tearing the form. Kits which have a torn AFTO FORM 104 will not be used to build the issuable kit until the MMM reinspects them.

Table 11. NSN 6545-01-400-3397 - Components First Aid Kit (Modified Kit of NSN 6545-01-094-8412)

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
6545-01-094-6142		CASE, MEDICAL INSTRUMENT AND SUPPLY SET, Nylon, Nonrigid, No. 8, 4-3/4 inches long x 2-3/8 inches wide x 5-1/2 inches high	1	EA
6510-00-159-4883		DRESSING, FIRST AID, Field, Camouflaged, 4 x 7 inches	2	EA
6510-00-200-3180		BANDAGE, GAUZE, Compressed and Folded, Brown, 2 inches x 6 yards	1	EA
6510-00-201-1755		BANDAGE, MUSLIN, Compressed, Brown, Camouflaged, 37 x 37 x 52 inches Triangular w/Safety Pins	1	EA
6530-01-316-0282		BOTTLE, Safety Cap, Plastic 16 Drams, 230s. Source of Supply: Kerr Group Inc., Distributor Don and Richards, 800-722-4866, Part Number CR-16, Cage OM7F6, \$33.00 CS of 230.	2/200	PG
		TOTAL WEIGHT TOTAL VOLUME		



Figure 7. First Aid Kit, Individual, NSN 6545-01-094-8412 (Sheet 1 of 2)

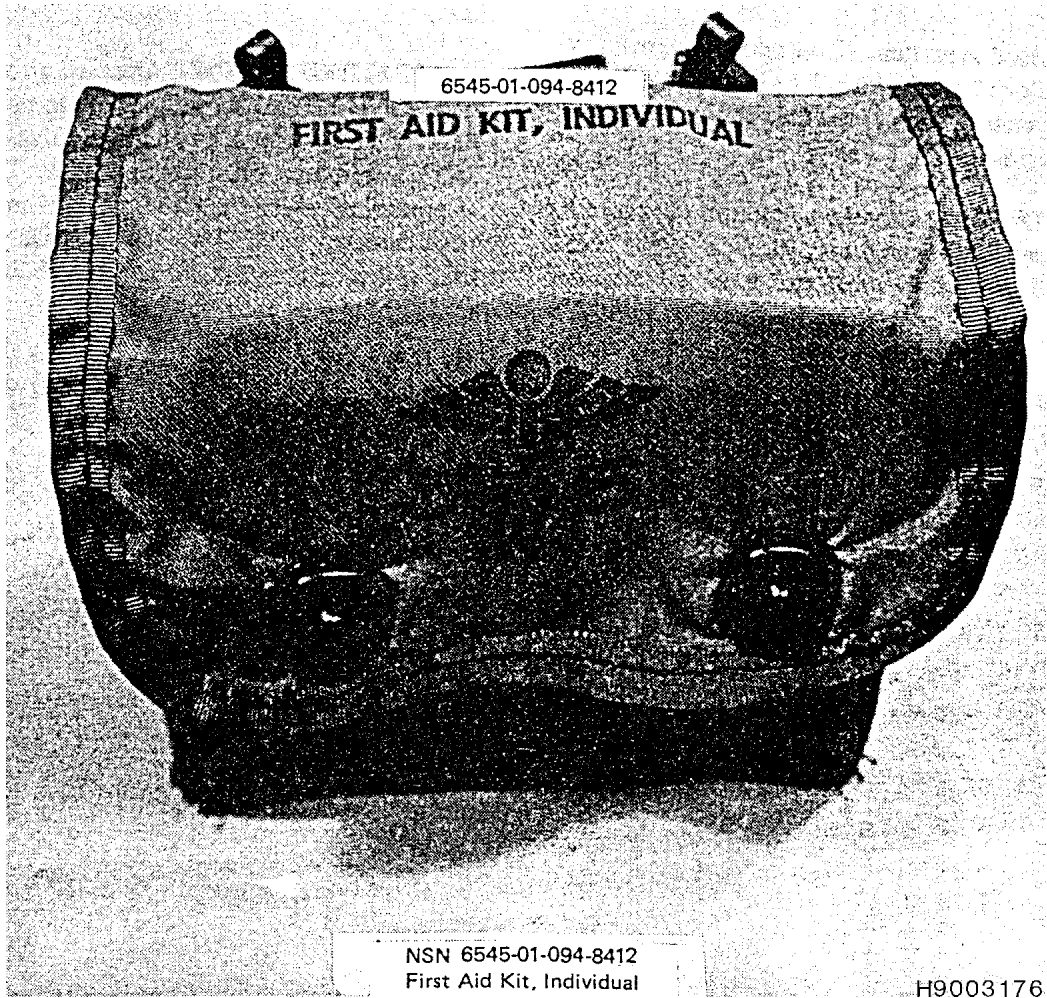


Figure 7. First Aid Kit, Individual, NSN 6545-01-094-8412 (Sheet 2 of 2)

23 NSN 6545-01-120-2632, SURVIVAL KIT, INDIVIDUAL, TROPICAL, TACTICAL AIRCREW, IN A PLASTIC BAG.

23.1 PURPOSE. This kit is specifically designed to provide a one week supply of first aid and survival supplies to stranded aircrew members during escape and evasion situations encountered under tropical conditions. This soft pack kit is mandatory for use with aircrew chemical defense equipment (ACDE) in accordance with T.O. 14S1-3-51. When ACDE is not required for use by aircrew members, this soft pack kit or NSN 6545-00-782-6412 can be installed in the SRU-21P survival vest, 12 pockets.

23.2 INSTALLATION. This kit will be issued in accordance with paragraph 14.2 of this TO.

23.3 INSPECTION. Using activities will return this kit to the MMM for inspection required by paragraph 7. The MMM will replace all "P"coded items that have less than two years of shelf life remaining unless items are not available. Follow instructions in paragraph 7.2 for "P"coded items with less than two years of shelf life remaining if they are needed to update the kit.

23.4 PACKING AND REPACKING.

23.4.1 Components of this kit are to be packed as flat as possible to keep the kit pliable for the survival vest. Aluminum foil strips are to be tucked under to prevent rupture of the inner plastic bag. Rubber bands are to be used to keep the strips folded. Items not encased in aluminum foil will be encased in a polyethylene bag as prescribed in paragraph 14.5.1. The completed kit packed in the first ziploc bag NSN 8105-00-837-7754 will be inserted (zippered end first) into the additional or second ziploc bag. This will provide more

durability for the kit while maintaining the size and flexibility of the current configuration. The completed/packed kit will not exceed 6 inches length.

Table 12. NSN 6545-01-120-2632 - Components

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
8105-00-837-7754		BAG, PLASTIC, Polyolefin, Interlock Seal Closure Fastening, 6 x 6 inches (Ziploc Bag), 1000s	1/1000	MX
4220-01-181-3154		FISHING KIT, EMERGENCY, Source of Supply: Fishing Kit and Insect Set Mittens Set can be purchased from Specialty Products Inc., 6626 Easton Rd., Pipersville, PA 18947, POC - Fred Prozzillo, Telephone 215-766-7777, CAGE OIRF7.	1	EA
6505-00-009-5060 (LP)	P	DOXYCYCLINE HYCLATE CAPSULES, USP, 100 mg, I.S., 100s, (Anti-Malarial Prophylactic Drug)	14/100	PG
6505-01-238-5632* (LP)	P	LOPERAMIDE HYDROCHLORIDE CAPSULES, 2.0 milligrams Loperamide, I.S., 100s (Anti-Diarrhea)	10/100	PG
6505-00-982-0288	P	ERYTHROMYCIN OPHTHALMIC OINTMENT, USP, 5 mg per gram, 1/8 oz (3.5 Gram)	1	TU
6510-00-054-7255		SKIN CLOSURE, ADHESIVE, SURGICAL, POROUS, 1/4 x 4 inches, 500s	10/500	BX
6510-00-913-7906		BANDAGE, GAUZE, ELASTIC, Two-Ply, 2 inches x 5 yards, 50s	1/50	BX
6510-00-913-7909	P	BANDAGE, ADHESIVE, 3/4 x 3 inches, 300s	12/300	BX
6510-01-483-3333	P	Towelette, benzalkonium chloride impregnated, 5" x 8", individually wrapped, 100s.	5/100	PG
6515-00-296-2529		SAW, Flexible Steel, 16 Teeth, 16 inches Long	1	EA
6515-01-363-1212		RAZOR, SURGICAL, PREPARATION, Small, Straight Type Single Blade, Plastic Folding Handle, 100s	1/100	PG
6695-00-118-2918 (LP)		BAG, FOOD SAMPLE, Polyethylene, Disposable, 18 fl oz Capacity, 500s	3/500	PG
6850-00-985-7166	P	WATER PURIFICATION TABLET, IODINE, 8 mg, 50s	1	BT
7610-01-301-5255 (LP)		LIST OF CONTENTS (Local Reproduction Authorized)	1	EA
*8415-01-192-2357		INSECT SET MITTENS SET, Head and Hand Covering, Nylon, Black	1	EA
8520-00-551-0375 (GSA Item)		SOAP, TOILET, Nonmedicated, Individually Wrapped, 1/2 oz, 1000s	1/1000	MX

Table 12. NSN 6545-01-120-2632 - Components - Continued

NATIONAL STOCK NUMBER	NOTES	DESCRIPTION	QTY	UNIT OF ISSUE
9920-01-191-3434 (LP)		MATCHES, Safety, Humidity Resistant, Waterproof Bag (1 Book of 20 Matches). Available from Opportunity Resources, Inc., POC Candace Lewis, 406-721-2930. Another possibility is NSN 9920-00-221-0613, Matches, Safety, 340s. Safety Matches can be purchased from GSA, or Datrex Inc., 3795 NW 25th St., Miami, FL 33142-6296, 305-638-8220. The Cage Code is 47823, Part Number 3561.	1	PG
		TOTAL WEIGHT 1.096 TOTAL VOLUME 0.041		

23.4.2 Items coded LP in the component listing may be procured from the manufacturer of the item or a local distributor. Before purchasing from a local source, check to make sure the item has not been recently depot stocked. Be aware that manufacturers and distributors have minimum order quantities.

23.5 MARKING. This kit does not require application of a red cross symbol.

23.6 SEALING AND/OR TAGGING. After completion of the inspection and repacking, the kit will be tagged with AFTO FORM 104 across the opening of the inner plastic bag. The AFTO FORM 104 will be affixed in a manner to preclude opening of the inner bag without tearing the form. Kits which have a torn AFTO FORM 104 will not be issued by life support or using activities. All such kits will be returned to the MMM for reinspection.

24 COMPONENT CROSS REFERENCE TO FIRST AID KITS.

This table is designed to cross reference a specific component to the kit NSNs it is in. The cross reference is in NSN component sequence.

Table 13. Component Cross Reference to First Aid Kits

COMPONENT NSN	KIT NSN	COMPONENT NSN	KIT NSN
H72678522	6545-00-116-1410	6510-00-200-5000	6545-01-090-0645
	6545-00-919-6650	6510-00-201-1755	6545-00-116-1410
	6545-00-922-1200		6545-00-919-6650
	6545-00-927-3000		6545-00-922-1200
	6545-01-090-0645		6545-00-927-3000
4220-01-181-3154	6545-00-139-3671		6545-01-090-0645
	6545-00-782-6412		6545-01-094-8412
	6545-01-120-2632		6545-01-400-3397
5335-00-373-2800	6545-00-116-1410	6510-00-201-3000	6545-00-116-1410
6505-00-009-5060	6545-00-782-6412	6510-00-201-7430	6545-00-927-3000
	6545-01-120-2632		6545-01-090-0645
6505-00-100-9985	6545-00-116-1410	6510-00-201-7680	6545-00-116-1410
6505-00-148-7096	6545-00-116-1410	6510-00-202-0800	6545-00-116-1410
	6545-00-922-1200	6510-00-203-5000	6545-00-927-3000
6505-00-254-5527	6545-00-927-3000	6510-00-203-5500	6545-01-090-0645
6505-00-491-7557	6545-00-116-1410	6510-00-597-7469	6545-01-090-0645
6505-00-926-2111	6545-00-116-1410	6510-00-782-2698	6545-01-090-0645
6505-00-982-0288	6545-00-139-3671	6510-00-913-7906	6545-00-782-6412
	6545-00-782-6412		6545-01-120-2632
	6545-00-853-6309	6510-00-913-7909	6545-00-116-1410
	6545-01-120-2632		6545-00-139-3671
6505-01-121-2335	6545-00-116-1410		6545-00-782-6412
6505-01-238-5632	6545-00-782-6412		6545-00-919-6650
	6545-01-120-2632		6545-00-922-1200
6508-01-436-0607	6545-00-116-1410		6545-01-094-8412
	6545-01-094-8412		6545-01-120-2632
6505-01-454-2363	6545-00-927-3000	6510-00-926-9082	6545-00-116-1410
	6545-01-090-0645	6510-00-935-5820	6545-00-922-1200
	6545-01-094-8412	6510-00-935-5822	6545-01-090-0645
6510-00-054-7255	6545-01-120-2632	6510-01-477-9736	6545-00-919-6650
6510-00-111-0708	6545-01-090-0645	6510-01-060-1639	6545-00-116-1410
6510-00-159-4883	6545-00-116-1410		6545-00-139-3671
	6545-00-919-6650		6545-00-782-6412
	6545-00-922-1200		6545-00-853-6309
	6545-00-927-3000		6545-00-919-6650
	6545-01-090-0645		6545-00-922-1200
	6545-01-094-8412	6510-01-107-7575	6545-00-853-6309
	6545-01-400-3397	6510-01-112-6414	6545-00-919-6650
6510-00-200-3075	6545-00-919-6650		6545-00-922-1200
	6545-00-922-1200	6510-01-483-333	6545-01-120-2632
	6545-01-090-0645	6510-00-296-2529	6545-00-139-3671
6510-00-200-3080	6545-00-116-1410		6545-00-782-6412
	6545-00-139-3671		6545-01-120-2632
6510-00-200-3180	6545-00-116-1410	6515-00-333-2200	
	6545-01-094-8412	6515-00-373-2100	
6510-00-200-3180	6545-01-400-3397	6515-00-383-0565	6545-00-927-3000
6510-00-200-3185	6545-00-919-6650		6545-01-090-0645
	6545-00-922-1200	6515-00-687-8052	6545-00-116-1410
6510-00-200-3190	6545-00-116-1410		6545-00-927-3000

Table 13. Component Cross Reference to First Aid Kits - Continued

COMPONENT NSN	KIT NSN	COMPONENT NSN	KIT NSN
6515-00-754-0426	6545-00-919-6650	6545-00-914-3500	6545-01-090-0645
	6545-00-922-1200	6545-01-094-6136	6545-01-094-8412
	6545-00-927-3000	6545-01-094-6142	6545-01-094-8412
6515-00-935-6592	6545-00-116-1410		6545-01-400-3397
6515-00-935-6593	6545-00-116-1410	6545-01-189-7094	6545-00-139-3671
	6545-01-090-0645	6545-01-189-7095	6545-00-782-6412
6515-00-935-7138	6545-00-116-1410	6605-01-326-1654	6545-00-782-6412
	6545-00-922-1200	6695-00-118-2918	6545-00-782-6412
	6545-01-090-0645	6840-00-142-8965	6545-00-782-6412
6515-00-958-2232	6545-00-116-1410	5955-00-879-7995	6545-00-116-1410
6515-01-225-4681	6545-01-090-0645	6850-00-985-7166	6545-00-116-1410
6515-01-363-1212	6545-00-782-6412		6545-00-139-3671
	6545-01-120-2632		6545-00-782-6412
6515-01-337-3689	6545-00-139-3671		6545-00-927-3000
	6545-00-782-6412		6545-01-094-8412
6520-00-926-9041	6545-00-116-1410		6545-01-120-2632
	6545-00-919-6650	7210-00-715-7985	6545-01-090-0645
	6545-00-922-1200	8105-01-195-8730	6545-00-927-3000
	6545-00-927-3000	8105-00-837-7754	6545-00-853-6309
	6545-01-090-0645		6545-01-120-2632
6530-01-316-0282	6545-00-116-1410	8315-00-787-7000	6545-00-927-3000
	6545-01-094-8412	8315-00-787-8000	6545-00-116-1410
	6545-01-400-3397		6545-01-090-0645
	6545-00-116-1410	8415-01-192-2357	6545-00-782-6412
6530-01-380-7309	6545-01-090-0645		6545-01-120-2632
6530-00-784-4205	6545-00-090-0645	8520-00-551-0375	6545-01-120-2632
6545-00-113-3722	6545-00-922-1200	8520-00-782-3554	6545-00-139-3671
6545-00-131-6992	6545-00-116-1410	8925-01-008-0960	6545-00-139-3671
	6545-01-120-2632	8935-00-082-5667	6545-00-139-3671
6545-00-853-6309	6545-00-116-1410	9535-01-201-7014	6545-00-139-3671
	6545-00-919-6650	9920-01-191-3434	6545-00-139-3671
	6545-00-922-1200		6545-00-782-6412
	6545-00-927-3000		6545-01-120-2632
	6545-01-094-8412		
6545-00-912-9830	6545-00-927-3000		
6545-00-912-9860	6545-00-919-6650		